

YOUR PARTNER IN PACKAGING



AC 1000

SEMI AUTOMATIC TABLET COUNTER



The AC 1000 offers small to medium scale counting operations with an unparalleled flexibility over a wide range of tablet types..



MACHINE FEATURES

- Counts coated and uncoated tablets, hard and soft gelatine capsules and dragees.
- Digitally controlled count adjustment.
- Accurate dispensing with excellent repeatability.
- Simple operation with fast product change over.
- Tool free removable product contact parts.
- Constructed in compliance with FDA approved materials.
- Optional full IQ/OQ Validation support documentation
- Small machine footprint

YOUR PARTNER IN PACKAGING

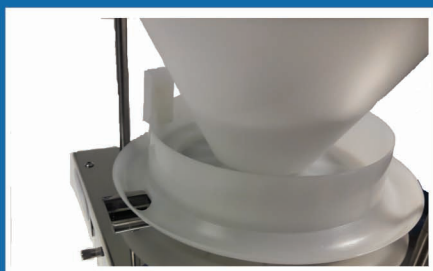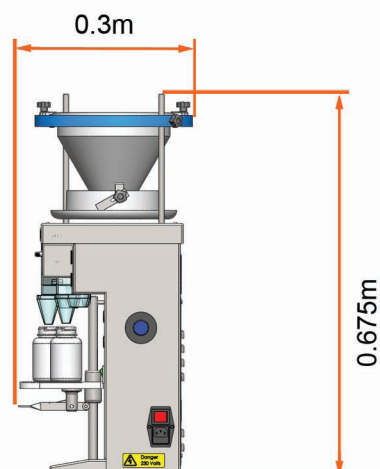

TECHNICAL SPECIFICATION

Output:	10 bottles per minute based on a 50 count Ø5mm round tablet
Electrical Control System:	Omron
Electrical Requirements:	220/240V 50Hz 1 Ph or 110/120V 60Hz
Build Standard:	cGMP, CE/UKCA marked in accordance with Machinery Directive 2006/42/EC
Product Size:	Tablet Ø5mm - Ø15mm Capsule #5 - #00
Container Size:	Diameter - 15mm 100mm Height - 25mm to 185mm
Nett Weight:	21Kg

OPERATION

Tablets or capsules are loaded into a product hopper which progressively feeds them into a vibratory scroll feeder to pre-orientate the product into a single stream.

The product proceeds to pass through a counting sensor where they are diverted through dispensing funnels into the pre-positioned container. When required count is reached the divider flap switches position and directs the product into the next container.